

[ SOLUTION SHEET: ]

<b>Solution:</b> Ascom Unite Axess
Smart device application for iOS and Android
Access to all messages, notifications and alarms when off site
Enables improved decision-making
High level of information security



## ASCOM UNITE AXESS

The Ascom Unite Axess is an enabling smart device application that helps healthcare providers and mobile clinicians deliver optimized patient care. It does so by effectively, securely and accurately handling and presenting alarms and other events, allowing caregivers access to prioritized and dedicated situation based information, anywhere and anytime.


To use Ascom Unite Axess, your healthcare facility must have previously purchased and installed Ascom Unite.

Within the modern healthcare industry, physicians, mobile clinicians and other caregivers increasingly want to use their personal smart devices at work. Ascom Unite Axess represents a powerful solution to this issue, enabling access to all messages, notifications and alarms when off site and functioning as a standard IP-DECT or VoWiFi handset on site.

Ascom Unite Axess enables complete integration with other Ascom solutions such as patient monitoring, nurse call, fire or other alarms and building management systems. Furthermore, the application also facilitates message diversion and escalation functionality, ensuring mission-critical information always receives an appropriate response, even if the first intended recipient is off site or does not respond.

Ascom Unite Axess also delivers high levels of patient and information privacy and security – all data is stored centrally, so even if a device is lost, data is safe. Furthermore, the application also offers support failover via SMS / SMTP to the same or another device and the increased security of end-to-end encrypted messages.

Operationally, Ascom Unite Axess delivers many benefits. Physicians and caregivers using the application are instantly more contactable. Quick, easy access to the right information or integration, regardless of location, enables improved decision-making and can have a positive effect on patient care optimization. Message prioritization presented via a smart, attractive and intuitive graphical user interface can also help physicians and other caregivers handle workloads more effectively.

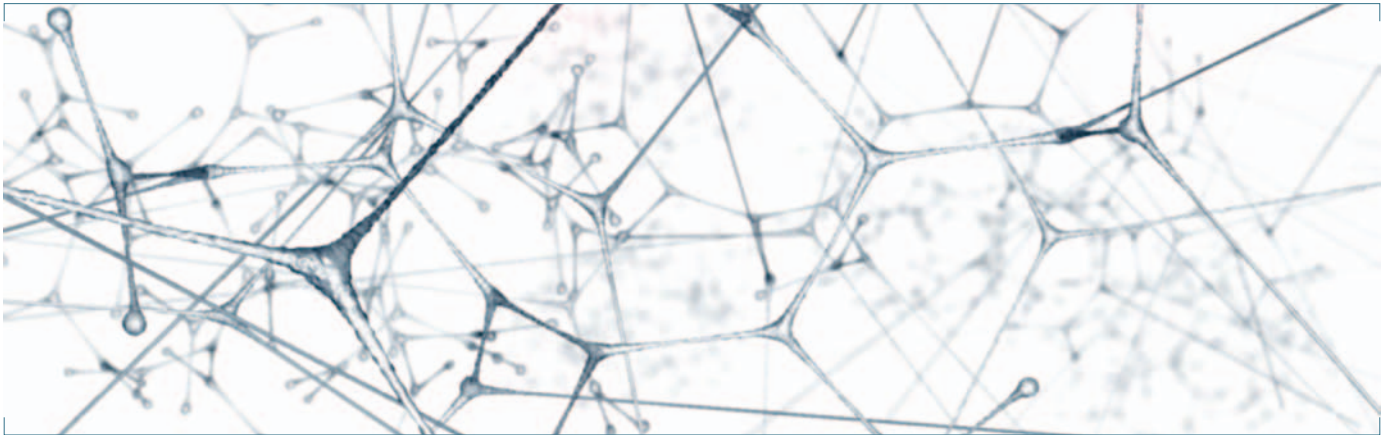

Enabling and secure, Ascom Unite Axess puts the potential to support improved standards of healthcare in the palm of your hand, wherever you are.

**Benefits Ascom Unite Axess:**

- Ascom Unite Axess is an enabling smart device application that helps healthcare providers optimize patient care.
- Enables complete integration with other Ascom solutions such as patient monitoring, nurse call, fire or other alarms and building management systems.
- Ascom Unite Axess enables complete Unite integration for messaging.
- Information is presented via a smart, attractive and intuitive Graphical User Interface.
- Diversions – if the user is not present on the WiFi or IP-DECT network, messages are automatically sent to the smart device.
- Escalations – if a message is not acknowledged within a pre-determined time, it is automatically escalated to another person, regardless of whether they are using VoWiFi, IP-DECT or a smart device.
- Intuitive message categorization offers easier access to certain types of message.
- Allows for the storage of all information on one device, independent of carrier technology.
- Centrally stored data prevents unauthorized data disclosure and optimizes security and patient integrity.
- Ascom Unite Axess produces an audit trail comprising of just one log for all messages and responses, independent of technology – this optimizes traceability and accountability.
- By using text messages and mail if no data connection is available, Ascom Unite Axess ensures notifications always reach the intended recipient, even if coverage is not optimal.